

shift³[®]

**Electronic Shiatsu
12-Point Foot
Massager**



USER MANUAL

IMPORTANT SAFETY INSTRUCTIONS

WHEN USING AN ELECTRICAL APPLIANCE, BASIC PRECAUTIONS SHOULD ALWAYS BE FOLLOWED, INCLUDING THE FOLLOWING:

READ ALL INSTRUCTIONS BEFORE USING THIS APPLIANCE.

DANGER - TO REDUCE THE RISK OF ELECTRIC SHOCK:


- ALWAYS unplug the appliance from the electrical outlet immediately after using and before cleaning.

WARNING - TO REDUCE THE RISK OF BURNS, FIRE, ELECTRIC SHOCK, OR INJURY TO PERSONS:

- An appliance should NEVER be left unattended when plugged in. Unplug from outlet when not in use, and before putting on or taking off parts.
- Do not operate under blanket or pillow. Excessive heating can occur and cause fire, electric shock, or injury to persons.
- Close supervision is necessary when this appliance is used by, on, or near children, invalids, or disabled persons.
- Use this appliance only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.
- Never operate this appliance if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water. Return the appliance to a service center for examination and repair.
- Do not carry this appliance by supply cord or use cord as a handle.
- Keep the cord away from heated surfaces.
- Never operate the appliance with the air openings blocked. Keep the air openings free of lint, hair, and the like.
- Never drop or insert any object into any opening.
- Do not use outdoors.
- Do not operate where aerosol (spray) products are being used or where oxygen is being administered.
- To disconnect, turn all controls to the off position, then remove plug from outlet.
- Use heated surfaces carefully. May cause serious burns. Do not use over insensitive skin areas or in the presence of poor circulation. The unattended use of heat by children or incapacitated persons may be dangerous.
- Do not stand on or in the appliance. Use only while seated.
- Do not use the product for more than 15 minutes at a time.

SAVE THESE INSTRUCTIONS

PRODUCT APPEARANCE and FUNCTION ILLUSTRATION

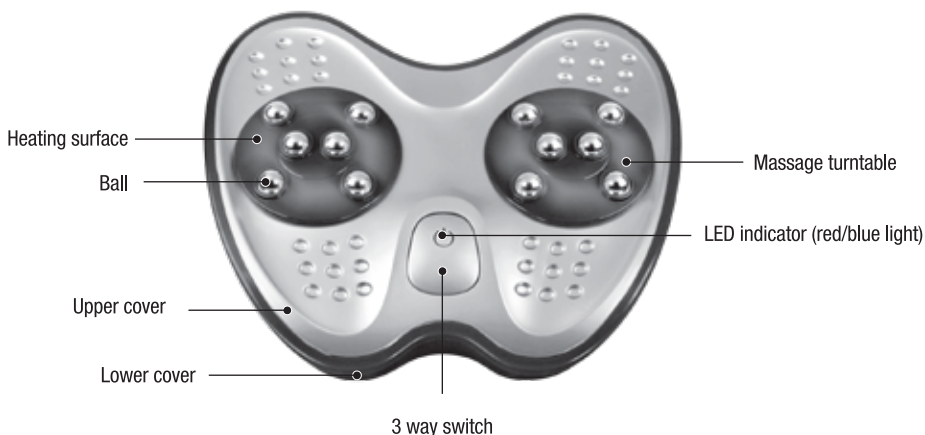

OPERATING INSTRUCTION

WARNING - To reduce the risk of electric shock, this appliance has a polarized plug (one blade is wider than the other). This plug will fit in a polarized outlet only one way. If the plug does not fit fully in the outlet, reverse the plug. If it still does not fit, contact a qualified electrician to install the proper outlet. Do not alter the plug in any way.

MAKE SURE THAT THE RATED VOLTAGE IS THE SAME OR COMPATIBLE WITH YOUR LOCAL ELECTRICAL VOLTAGE.

3 way switch: massage, massage with heat, OFF.

- (Kneading Massage): Press the 3 way switch button for the first time. The LED light will be BLUE. The two turntables will start rotating in mutual reverse directions.
- (Kneading and Heat Massage): Press the 3 way switch button again. The LED light will be RED. The two turntables will start rotating in mutual reverse directions. The lights on the turntables will be lit up and the surface of the turntables will generate heat.
- Power OFF: Press the 3 way switch button for the third time. The LED light will be off, the two turntables stop rotating and the lights on them extinguish as well.

CLEANING

- Unplug the product and coil the cord into the storage space on the bottom. Place the product on a flat, stable surface. Use a soft and damp sponge only. Never drop water or any other liquid onto the massager.
- Do not immerse the product into water for cleaning.
- Never use abrasive cleaner, brush, gasoline, kerosene, furniture glazing, coating or thinner to clean the product.

STORAGE

- Put the product into package when not in use and store neatly in a dry, cool location protected from direct sunlight.
- Avoid sharp edges or pointed objects, as they may puncture the surface of the product.
- Do not use the power cord to carry the massager.

MAINTENANCE

- If a failure or fault occurs, the product should be returned to either the retailer from where the item was purchased for replacement or return. Do not repair by yourself.

TECHNICAL SPECIFICATION

- Power Supply: AC 120V/60Hz.
- Power Consumption: 25W

Foot Massager MerchSource, LLC	CONFORMS TO UL STD. 1647 CERTIFIED TO CAN/CSA STD. C22.2 No. 68-92
	For Household Use Only. Pour usage domestique seulement Use Indoors Only Utiliser seulement à l'intérieur.
MODEL 1618395 120V, 60Hz, 25W	3173063 Made in China/Fabriqué en Chine

"SHIFT 3" and the SHIFT 3 LOGOS and other indication of origin are registered trademarks of MerchSource, LLC. All logos, text and graphics are also the subject of copyright ©2009 of MerchSource, LLC. All rights reserved.

©2009 MerchSource, LLC. All Rights Reserved.
Distributed by MerchSource, LLC, Foothill Ranch, CA 92610.

The artwork and design of this manual are protected by US copyright law and may not be reproduced, distributed, displayed, published or used for any purpose without prior written permission. It is not permitted to alter, remove or reproduce any trademarks or copyright notice from this package.